



UV-B Narrowband PL-L/PL-S

PL-L 36W/01/4P 1CT/25

More than 400 independent clinical studies have proven that the UVB Narrowband treatment is safer and more effective than any other treatment in its class. Lamps installed in such phototherapy treatment systems emit only a very narrow waveband from the 'B' bandwidth of the UV spectrum (290 to 315). Philips offers lamps with narrow waveband of between 305 and 315 nm which peaks at 311 nm. This makes these lamps very suitable for UV-B Narrowband phototherapy systems which treat skin diseases such as psoriasis and vitiligo. The PL-L/PL-S versions provide additional design freedom since these are space-saving compact, single-ended lamps. Further flexibility is assured since these lamps use the same lamp caps as general lighting lamps and also use the same universal ballasts.N.B.: Our UVB lamps are NOT registered with FDA as medical devices as they are NOT packaged or labeled for commercial distribution for health-related purposes. US customers are referred to the UVB and UVA lamp range brochure US version.

Product data

General Information			
Cap-Base	2G11 [2G11]		
Main Application	Phototherapy Systems		
Life to 50% Failures Preheat (Nom)	1000 h		
Useful Life (Nom)	1000 h		
Light Technical			
Color Code	01		
Chromaticity Coordinate X (Nom)	228		
Chromaticity Coordinate Y (Nom)	230		
UV Depreciation at 500 h	15 %		

UV Depreciation at 1000 h	20 %
Operating and Electrical	
Power (Rated) (Nom)	36 W
Lamp Current (Nom)	0.435 A
Voltage (Nom)	105 V
Mechanical and Housing	
Cap-Base Information	4 Pins
UV	
UV-B Radiation 100 hr (IEC)	4.8 W

Datasheet, 2018, September 14 data subject to change

UV-B Narrowband PL-L/PL-S


UV-B Radiation 5hr (IEC)	6.2 W
Product Data	
Full product code	871150086889340
Order product name	PL-L 36W/01/4P 1CT/25
EAN/UPC - Product	8711500868893
Order code	927903400121

Numerator - Quantity Per Pack	1	
Numerator - Packs per outer box	25	
Material Nr. (12NC)	927903400121	
Net Weight (Piece)	104.000 g	

Warnings and Safety

• A lamp breaking is extremely unlikely to have any impact on your health. If a lamp breaks, ventilate the room for 30 minutes and remove the parts, preferably with gloves. Put them in a sealed plastic bag and take it to your local waste facilities for recycling. Do not use a vacuum cleaner.


Dimensional drawing

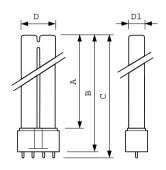


PL-L 36W/01/4P

Product	D1 (max)	D (max)	A (max)	B (max)	C (max)
PL-L 36W/01/4P 1CT/25	18 mm	39 mm	384.2 mm	410 mm	416.6 mm

Photometric data

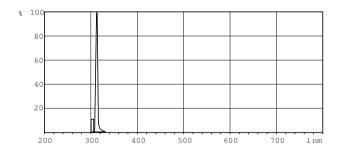

Lightcolor /01

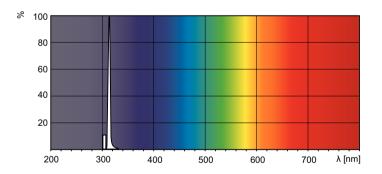

Lightcolor /01

UV-B Narrowband PL-L/PL-S



© 2018 Signify Holding All rights reserved. Signify does not give any representation or warranty as to the accuracy or completeness of the information included herein and shall not be liable for any action in reliance thereon. The information presented in this document is not intended as any commercial offer and does not form part of any quotation or contract, unless otherwise agreed by Signify. Philips and the Philips Shield Emblem are registered trademarks of Koninklijke Philips N.V.